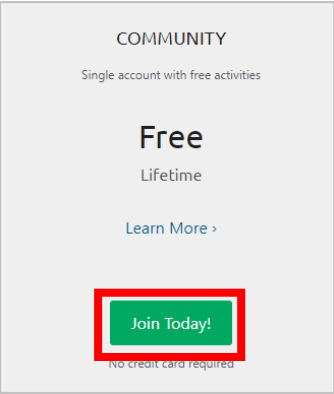
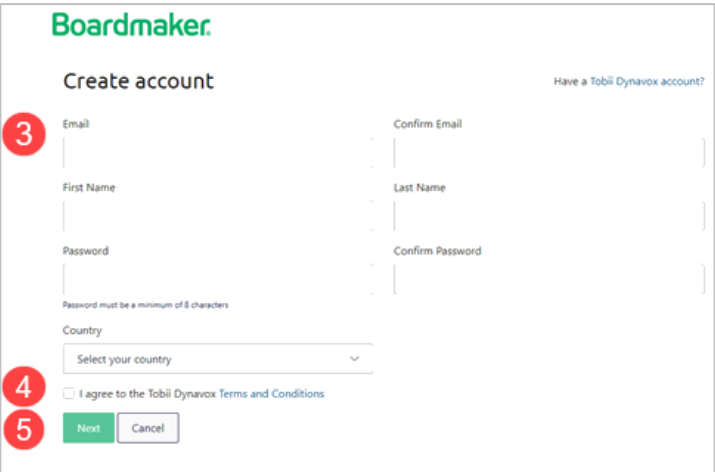


Create a Free Boardmaker Community Account: Quick Guide

A Boardmaker Community account allows you to search for Boardmaker activities created by other Community members and the Boardmaker Education Team, download activities to edit in your Boardmaker software, join groups specific to your interest, and share Boardmaker Activities you create. If you use other Boardmaker or Tobii Dynavox products, the Community is even more powerful. Sync your Community account to your Boardmaker 7 Standard license to store activities in the cloud. Add Community activities to the Snap Activities folder and play them through the Boardmaker Player in TD Snap.

This guide provides step by step instructions for creating and getting started with a free Boardmaker Community account.

Create a Boardmaker Community Account

1. Go to get.myboardmaker.com. 2. In the Community column select Join Today!	
3. Fill out the form. 4. Read the terms and check the box if you agree 5. Select Next.	

6. Read the terms and check the box if you agree.
7. Select **Start**.

Get Started in the Boardmaker Community

1. Complete your Community profile. Edit your username, organization, occupation, and add a picture if desired.



NOTE: Edit your profile any time by selecting your name in the top right corner → My Account.

2. Select Standards.
 - a. Select the plus sign to expand the standards.
 - b. Check the box to add the standards.



NOTE: You can search Community activities by the standards you select.

3. Select **Next**.

Learn More

Learn more about the Boardmaker Community and other Tobii Dynavox products at learn.tobiidynavox.com.

The Boardmaker Player in TD Snap – Adding Activities

The integrated Boardmaker Player in TD Snap allows you to play preloaded and interactive activities from Boardmaker right in your TD Snap software. Preloaded activities include an accessible bookshelf, accessible games, an interactive calendar, math tools, as well as samples from Boardmaker's Core First Learning, Reading Avenue, and Expedition Education curriculum titles. This document will explain how to add activities from your Boardmaker account to be used in the integrated Boardmaker Player in TD Snap. Select the statement below that best describes your need to access the appropriate steps to add activities to the TD Snap Boardmaker player.



myTobiiDynavox
+
Boardmaker

1. The device user's TD Snap is signed into myTobiiDynavox using the same email address that is associated with the paid Boardmaker account that will be sharing activities.

Example: a parent has a paid Boardmaker account and uses the same email address for their communicator's myTobiiDynavox account. If they are signed into myTobiiDynavox in the TD Snap software, they can share activities to the communicator's integrated Boardmaker Player following [these steps](#).



myTobiiDynavox



Boardmaker

2. The email address associated with the paid Boardmaker account sharing activities is different than the myTobiiDynavox account signed into the communicator's TD Snap software.

Example: a teacher or therapist has a paid Boardmaker account and supports multiple communicators using TD Snap. In this scenario, the teacher or therapist will need to share the activity to the communicator's Boardmaker account first. To share activities




Scenario 1

The email address for myTobiiDynavox and the Boardmaker paid account are the same.

1. Log into your Boardmaker account to create or locate activities to add

- a. Visit myboardmaker.com
or
- b. If using a device with internet, navigate to the Dashboard in TD Snap and access the **Boardmaker Library**
- c. Select the **Go to the Boardmaker** Icon 
- d. Select **Log In**
- e. Log in using your Tobii Dynavox single sign on email address and password

Note: Use the same email address that is logged into myTobiiDynavox in TD Snap in order to connect the two accounts
- f. Select **Log In**

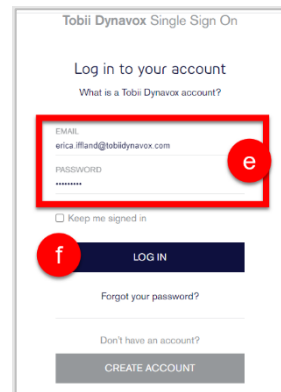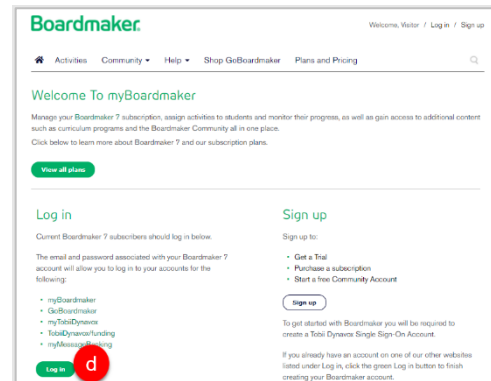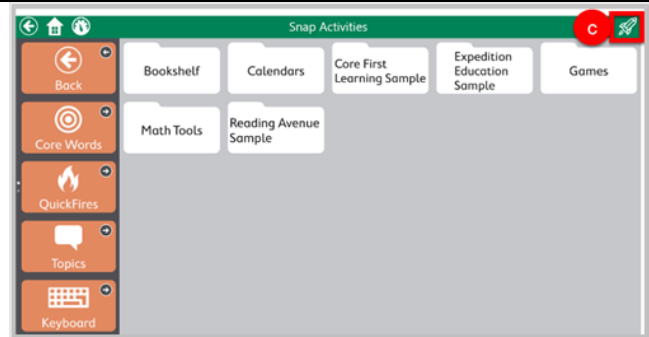


2. Add activities from the Community to My Boardmaker

Note: If you created activities in the Boardmaker Editor, or opened them through the email functions in myBoardmaker.com, skip this step and go to step 3

Search the Community


- Click on the **magnifying glass** in the menu bar
 - Type in your search term (i.e.: activity type, topic, etc.)
 - Hit enter on your keyboard
 - Scroll through or use filters to find activities
 - Select **the Title** of the activity to see a preview
 - Select the **+** to add the activity to My Boardmaker
- Note:** The blue tool bar will change when the activity is successfully added to My Boardmaker
- Once successfully added, select the **My Boardmaker** tab

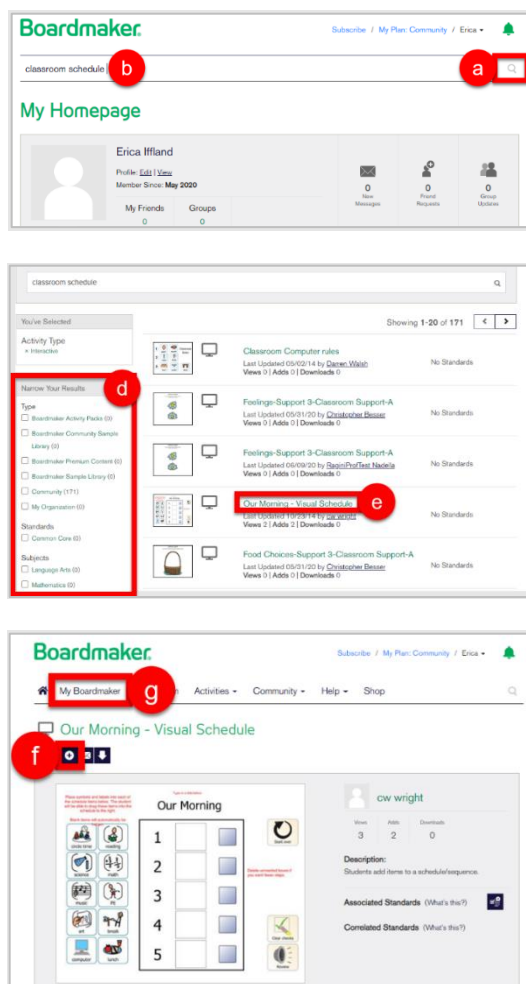


3. Add activities from My Boardmaker to the integrated Boardmaker Player


Add any activity to the Snap Activities Folder

- From the My Boardmaker tab, locate the activity you want to add
- Select the **Move activity to folder** icon 

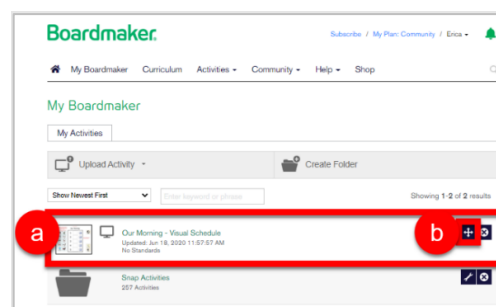

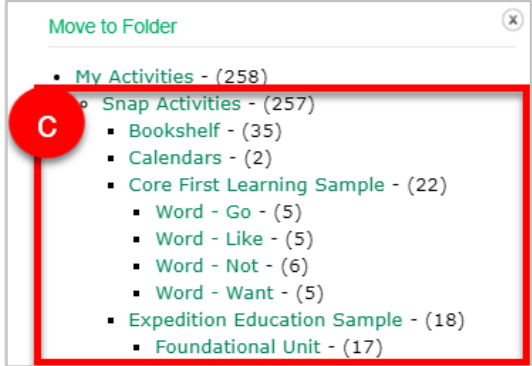c. From the popup, select Snap Activities or an appropriate subfolder within Snap Activities	
--	--

4. Playing Activities in the Boardmaker Player

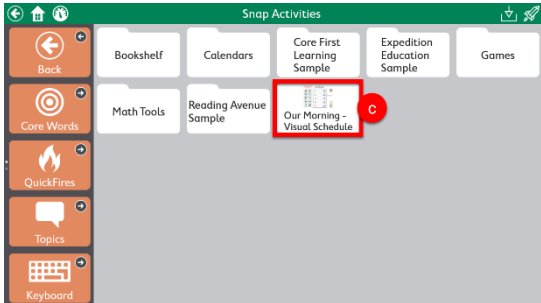a. Open TD Snap b. Navigate to the Dashboard and select Boardmaker Library c. If connected to Wi-Fi and logged into the associated myTobiiDynavox account in TD Snap, your activity is now available to be used in the integrated Boardmaker Player Note: If the Boardmaker Library was already open in TD Snap you may need to navigate out and reopen to see the addition of the activity	
--	---

Note: Any activity added to the Snap Activity folder in My Boardmaker will appear in the integrated Boardmaker Player of any device that is logged into the myTobiiDynavox account using the same email address. If you are working with multiple AAC users or want to assign specific activities to each user, see the next section.

Scenario 2

The email address for myTobiiDynavox and the Boardmaker paid account are different.

Before you begin:




- The individual (i.e., teacher or therapist) creating the activities for the TD Snap user must have a Boardmaker subscription that allows them to create interactive Boardmaker activities.
- The TD Snap user must create a free Boardmaker community account. Learn how by using this [guide](#).

Note: Make sure you create the free Boardmaker account using the same email and password you use to log into your myTobiiDynavox account within the TD Snap software

- The TD Snap user must share the email address they used to set up their account with their teacher or therapist who will be using their Boardmaker account to create and share activities.

1. The teacher or therapist with a paid Boardmaker account shares activities with the TD Snap User

- Visit <https://www.myboardmaker.com/login.aspx>
- Select **Log In** and log into your (teacher/therapist) paid account
- Select the **My Boardmaker** tab
- Select the activity you want to add to the communicator's Boardmaker Player in TD Snap
- Select the **Mail** icon 
- Type in the communicator's email address (the email address used to create their myTobiiDynavox and Boardmaker Community accounts)
- Select **Share**
- Repeat this process if sharing multiple activities
- Log out of your account


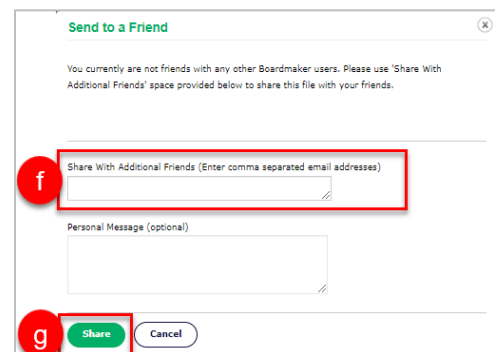
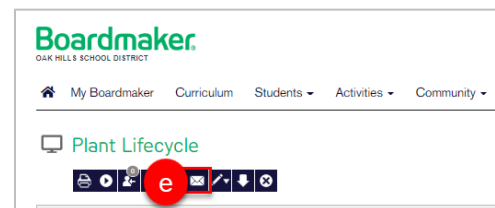
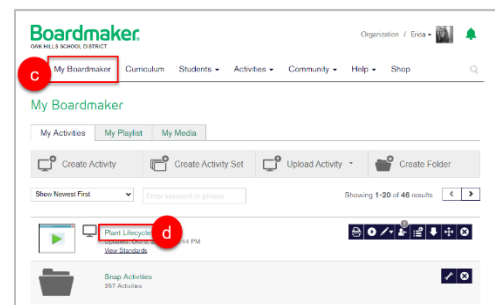


Before you begin step 2:

- You must have the log in information for the TD Snap user's Boardmaker Community account to complete this step
- Check to be sure the student is signed into myTobiiDynavox in their TD Snap software (you will need internet access for this step)
 1. Open **TD Snap**
 2. Select the **edit icon** in the top right corner
 3. Select **User**
 4. Select **Log In** under the myTD Account heading
 5. Either input the email and password or select Create Account
 6. Follow prompts on the screen
 7. Once logged in select **Done** to exit edit mode in TD Snap

2. The teacher, therapist, or family adds the shared activities to the communicator's Boardmaker Player in TD Snap



- Visit <https://www.myboardmaker.com/login.aspx>
 - Select **Log In** and log into the TD Snap user's Boardmaker Community account
 - Select **New Messages** from the home screen
 - Select the message Subject "**Check Out A File From a Friend...**"
 - Select the **Activity Link**
 - Select the **+ button** to add to your myBoardmaker
- Note:** return to your home screen and repeat steps c - f if multiple activities were shared
- Select the **My Boardmaker** tab
 - Find the activity and select the **move activity to folder icon** 

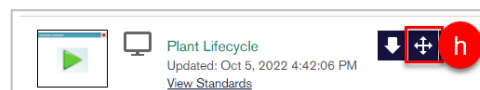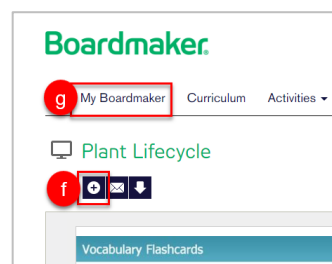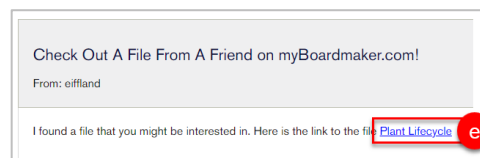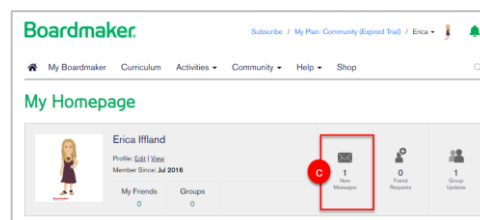

- i. Select the **Snap Activities** folder

Note: repeat steps h & i if multiple activities were shared

- j. Open **TD Snap** on the communicator's device
- k. Select the **Dashboard**
- l. Select **Boardmaker Activities**
- m. Your activity will now be available within the Boardmaker Library
- n. Play the activity and have fun!

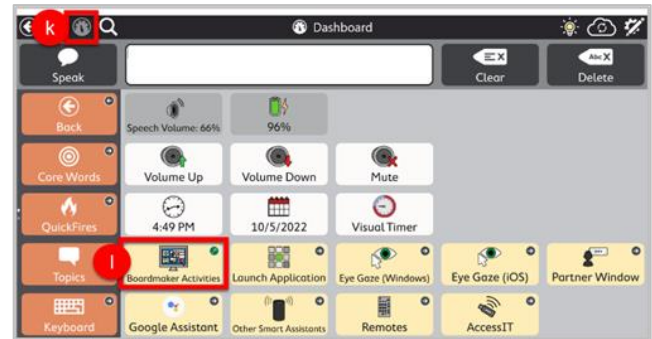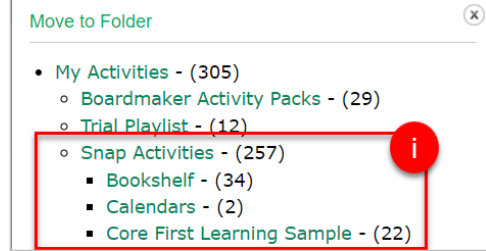


Explore Pageset Central for TD Snap®

Pageset Central is a collaborative community platform on www.mytobiidynavox.com where AAC users and professionals can share and download Page Sets and Page Bundles for Tobii Dynavox communication software — supporting the unique needs and preferences of AAC users beyond the standard tools Tobii Dynavox offers.

What Types of Resources are Available on Pageset Central?

Within Pageset Central, you will encounter two primary types of resources: Page Sets and Page Bundles. While both are integral to customizing your communication experience, they serve different purposes and are uniquely structured.


- **Page Set:** A Page Set is a complete, standalone layout ready for immediate use upon upload to TD Snap. It includes a carefully arranged set of buttons, actions, and navigation tools crafted to support communication effectively.
- **Page Bundle:** A Page Bundle either contains a singular page or a group of pages. These bundles are uploaded to TD Snap and integrated within a Page Set, effectively acting as a toolkit for various environments, interests, and communication needs.

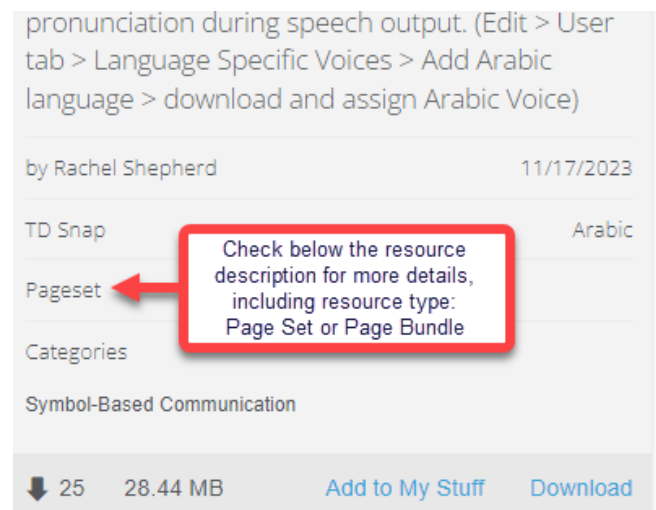

Note: Although Page Sets and Page Bundles are uniquely formatted, the steps to search, download, and import them are nearly identical. These terms will be used interchangeably in the remainder of the Quick Guide, unless otherwise specified.

Getting Started on Pageset Central

1. Confirm you are connected to the internet.
2. Navigate to www.mytobiidynavox.com in your internet browser.
3. If you have created an account with Tobii Dynavox or Boardmaker on one of our other websites or software, log in with the same account (SSO) or sign up for a free Tobii Dynavox account.
4. Select Pageset Central.

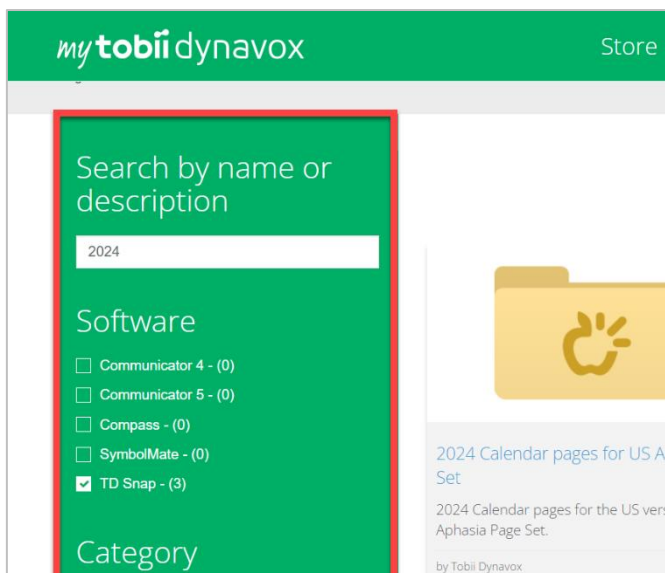

Searching for Page Sets and Bundles

Use the search bar on the left-hand side of the screen to search Page Sets by name or description.


You can also use the filters to narrow down search results. Available filters include:

- **Software** – filter by your preferred TD communication software.
- **Category** – filter by purpose, intended use, or content.
- **Language** – filter by preferred language.

Downloading Page Sets and Bundles

Once you've located the Page Set you want to download from Pageset Central, there are two ways you can download it:

1. If you are logged into your mytobiidynavox account, you can select "Add to My Stuff". This will download the Page Set to your mytobiidynavox account, making it available to import into TD Snap later via your synced account.
2. If you prefer to download the Page Set to a USB or directly to your computer files instead, you can select "Download." Once downloaded, you can either move the file onto a USB or keep it in your computer files – this will allow you to import into TD Snap® later.

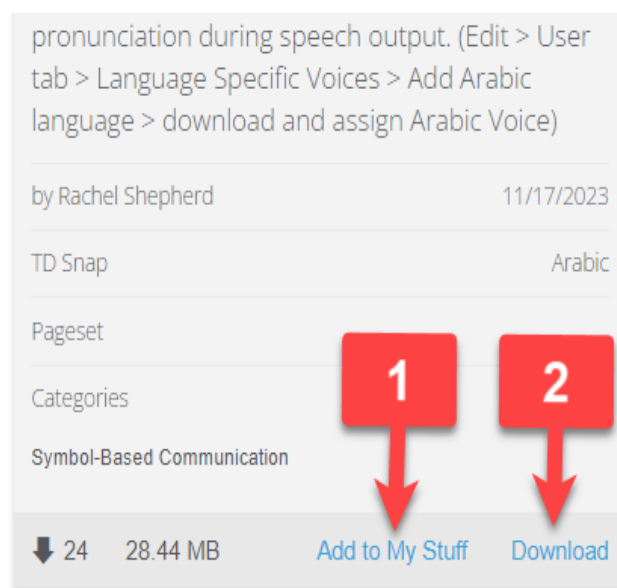

Importing Page Sets and Bundles into TD Snap®

Open TD Snap on your device and confirm you are connected to the internet before proceeding through the steps below.

Importing a Page Bundle

1. Confirm you are in the Page Set you want to add the Page Bundle to (e.g., Aphasia Page Set, Core First, Motor Plan).
2. Next, open the page from which you want the new Page Bundle to link to (e.g., Keyboard page).
 - **Tip:** For example, if you want to link the Emoji Page Bundle to a button in the Keyboard page in TD Snap Motor Plan, you will want to open the Keyboard page before moving on to Step 3.

3. **Select the Edit icon (gear and pencil)** in the upper right corner of the app.
4. Select the button or cell that you would like to link the Page Bundle to. You may have to scroll down to find an empty button or cell.
5. Select the **Button** tab.
6. Select **Link to Page** from the left-hand column.
7. Select **Import and Link to Page**. Depending on where you downloaded the Page Bundle to, you will select from either:
 - mytobiidynavox (you'll be prompted to log in to your account if you aren't already)
 - or a local Page Bundle File (you'll be prompted to browse your USB or computer for a .spb file)
8. Select the Page Bundle you want to import then **Next**.
9. Select the page you want this button to link to and then select **Import**.
 - **Tip:** Some page bundles may contain Topic or Main List pages (e.g., Main List: Food, Topic: Fast Food) with subordinate pages (e.g., Related List: Food – Breakfast). If prompted with this option, **you will typically want to pick the Topic or Main List**.
10. Select **Done** in the upper right-hand corner to try it out.

Importing a Page Set

1. With TD Snap open, **select the Edit icon (gear and pencil)** in the upper right corner of the app.
2. Select the **Page Set** tab.
3. Select the name of the current Page Set (e.g., Aphasia Page Set, Core First, Motor Plan).
4. Depending on where you downloaded the Page Set to, you will **either select Create Page Set or Download Page Set**.
 - If you downloaded the Page Set to your **Computer or USB**, select **Create Page Set** and continue with the instructions below to finish importing your Page Set into TD Snap.
 - a. After you select **Create Page Set**, select **Local Files**.
 - b. Browse your local files for a Page Set file (.sps).
 - c. Select **Next**.
 - d. You will have the option to rename the Page Set – this is optional.
 - e. Select **Create**.
 - f. Select the newly created Page Set from your list in the Page Set tab.
 - g. Select **Done** in the upper right-hand corner to try it out.
 - If you added the Page Set to **“My Stuff” or your mytobiidynavox account**, select **Download Page Set** and continue with the instructions below to finish importing your Page Set into TD Snap®.
 - a. After you select Download Page Set, select the Page Set you want to download from the list and then **select Download**.
 - b. You will have the option to rename the Page Set – this is optional.
 - c. Select **Save**.
 - d. Select the Page Set you downloaded from the list of Page Sets.
 - e. Select **Done** in the upper right-hand corner to try it out.


Uploading Page Bundles from TD Snap® to Pageset Central

Are you interested in uploading a Page Bundle you created in TD Snap to Pageset Central? Start by naming each TD Snap page you want to share with clear and easily searchable terms so others can easily find your Page Bundle on Pageset Central.

For example, if you're creating a book, consider naming each page with a keyword from the book title along with the page number. For example: PeteCat 1, PeteCat 2, and so on.

You can rename pages in TD Snap by:

1. Enter edit mode. 
2. Select the **Page** tab.
3. Type the desired page title into the **Label** field.

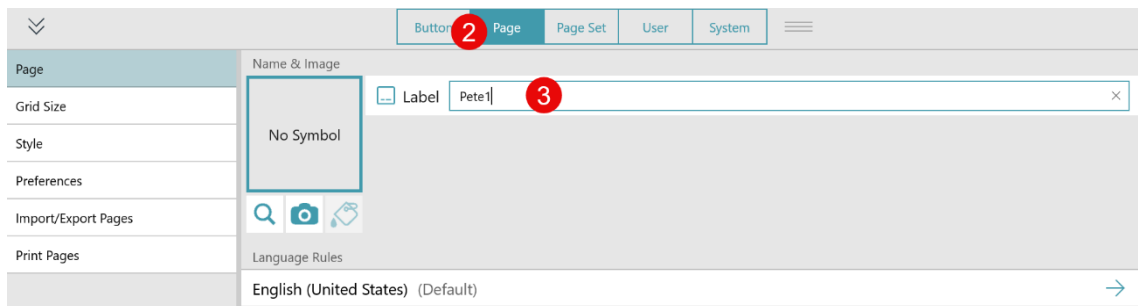

4. Select **Done** to exit edit mode.




Once you have the page titles updated, confirm you are connected to the internet and logged into your mytobiidynavox.com account in TD Snap. Then, proceed through the steps below.

In TD Snap®


1. Enter edit mode. 
2. Select the **Page** tab.
3. Select **Import/Export Pages**.
4. Select **Export Page Bundle**.
5. Select **To myTobiiDynavox**.



6. **Search** for the title of the first page in your bundle.
7. **Select** all additional pages included in the bundle (e.g., all the pages that can be navigated to from you're the first page you created).
8. When all necessary pages are selected, choose **Next**.

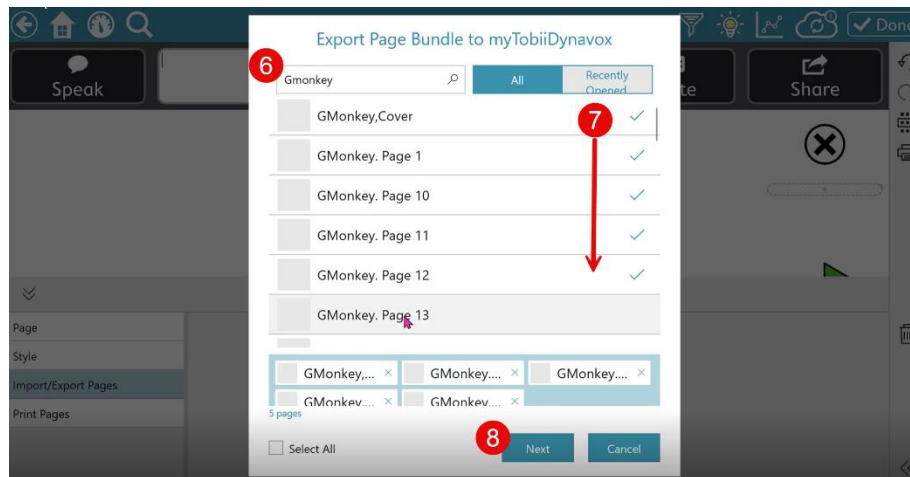


9. Name your Page Bundle and provide a description. Make it clear enough that others will know exactly what they will be importing if they choose to download it.
10. Select **Export**.
11. When it's completed, select **Done**. It is now saved in your personal mytobiidynavox.com account.

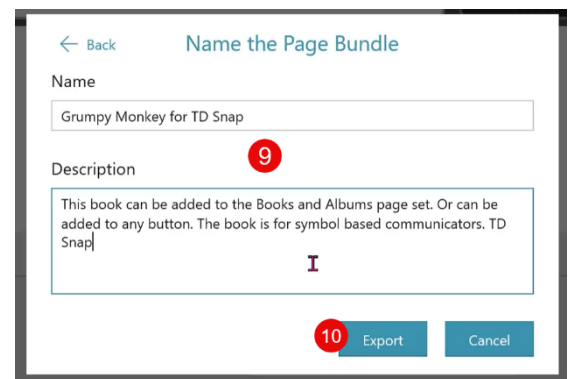


In myTobiiDynavox.com


1. Log in to your account.
2. Select the TD Snap tile.
3. Select the **Pagesets** tab.
4. Find the Page Bundle you uploaded. Select **Manage**.

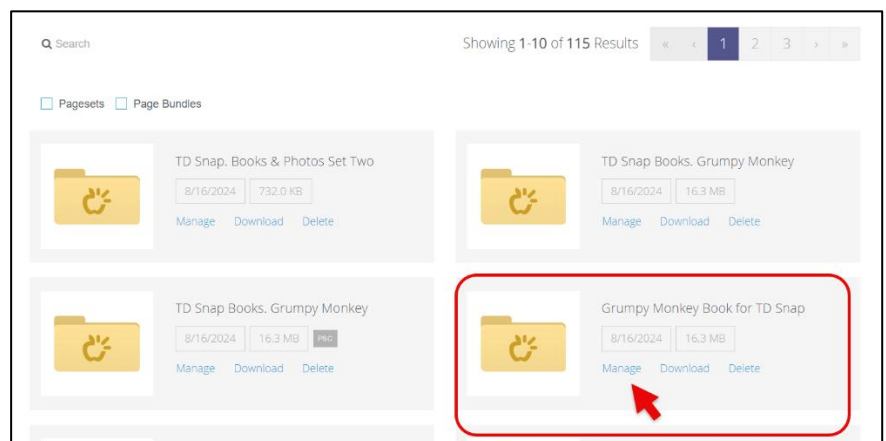


5. Select **Edit/Share to PSC**.

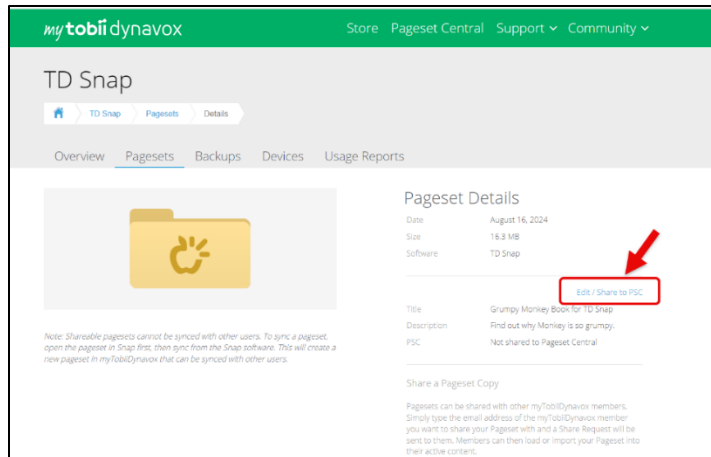

6. Toggle on **Share to Pageset Central**.

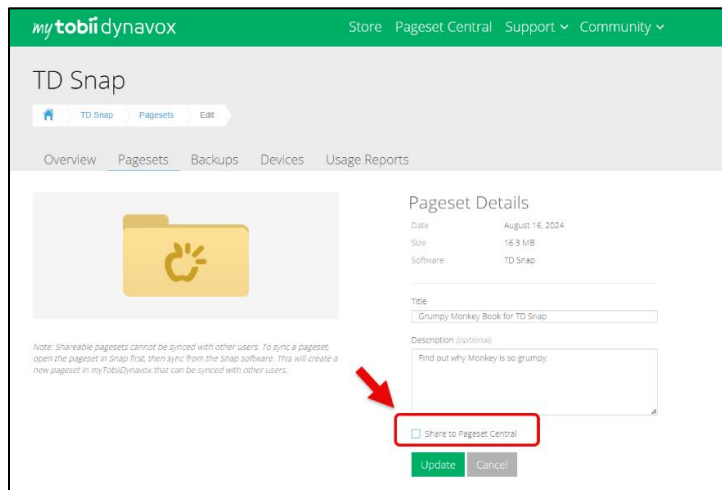

7. Confirm your title, description, and add a name you'd like to publicly display as author. Select **Update**.

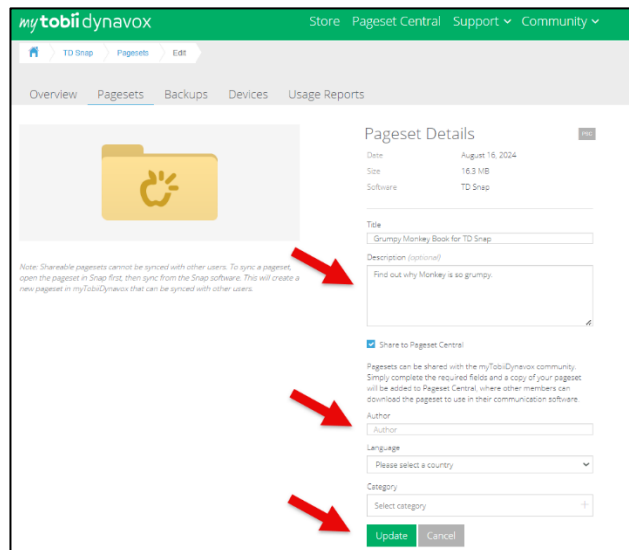

Continue Exploring Pageset Central

Looking for examples to kickstart your search on Pageset Central? Explore a selection of user-created Page Sets and Page Bundles below, providing a glimpse into the diverse options available for download.

Page Sets

- [TD Snap Core First English - Chinese Bilingual Page Set](#) and other Page Sets designed for people who need multilingual AAC supports.
- [Hospital Rehab Page Set](#) — adapted from the Aphasia Page Set and designed for people who need AAC support in a hospital setting.

Page Bundles

- [Rating and Pain Scales](#) Page Bundle
- [Emojis](#) Page Bundle — formatted for use in TD Snap® Motor Plan 66 on devices using the iOS platform.
- [Peppa Pig](#), [Minecraft](#), [Disney and Barbie](#), and other themed Page Bundles.

Continue Your Learning

- Access the [TD Snap® User Manual](#) to learn more about TD Snap®.
- Explore our [live calendar](#) for more virtual trainings on topics you care about led by our experts.
- Learn anytime by logging into the [Tobii Dynavox Learning Hub](#) for free resources like videos, guides, implementation tips, and more.
- Contact us by reaching out to your local [Solutions Consultant](#) for support with hardware, software, device funding, and more.